



MODIFICATION UNITS

R2110

Variable Speed, Single Spindle Unit

DESCRIPTION: Electric modification unit (custom built by Art Optical) with variable speed function from 0-1550 rpm's. The motor is housed in a Formica covered box. The tapered shaft extends into a brushed stainless steel pan that is made water-tight through the use of rubber o-rings. The pan is removable for cleaning.

INCLUDES: R2300 - Cone shaped sponge with plexiglass base and screw set
 R2260 - Suction cup (2 piece)
 R1550 - Sil-02 Care Polish (4 oz.)
 R2307 - Round 3"x1" replacement sponge with plexiglass base and screw set
 Allen key to remove tapered arbor
 Instructions for modification unit

R1590

Boston Manual Polishing Kit

DESCRIPTION: Hand held circular hoop (like stitching hoop) that has velveteen stretched across it. Modifications are done by hand, using Boston Professional Cleaning Polish and a two-piece green suction cup.


INCLUDES: Manual polishing hoop
 R2700VP - Velveteen pad (5")
 R1500 - Boston Professional Cleaning Polish (2 oz.)
 R2260 - Suction cup (2 piece)

R2100

Single Speed, Single Spindle Unit

DESCRIPTION: Same as variable speed unit except there is one constant spindle speed of 1550 rpm's.

INCLUDES: R2300 - Cone shaped sponge with plexiglass base and screw set
 R2260 - Suction cup (2 piece)
 R1550 - Sil-02 Care Polish (4 oz.)
 R2307 - Round 3"x1" replacement sponge with plexiglass base and screw set
 Allen key to remove tapered arbor
 Instructions for modification unit

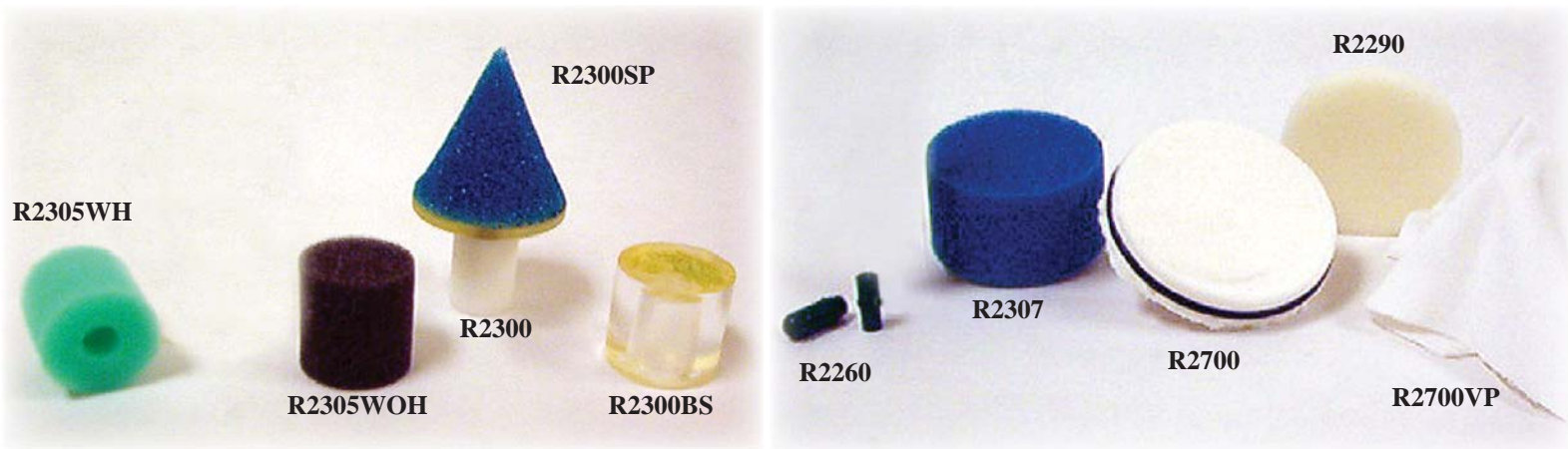

MODIFICATION ACCESSORIES

R2305WH

Sponge with Hole (Replacement)

DESCRIPTION: Replacement sponge for polishing. The plastic polishing tool base does not have to be returned. The old sponge can be removed from the tool base and the replacement sponge inserted.

R2305WOH

Sponge without Hole (Replacement)

DESCRIPTION: Replacement sponge for polishing. The plastic polishing tool base does not have to be returned. The old sponge can be removed from tool base and the replacement sponge inserted.

R2300

Cone Shaped Sponge with Polishing Tool Base

DESCRIPTION: Cone shaped sponge screw mounted on a plastic polishing cone tool base. Used to clean and polish the inside (base curve side) of gas permeable and PMMA lenses.

R2300SP

Cone Shaped Sponge (Replacement)

DESCRIPTION: Replacement cone shaped sponge for polishing cone tool base.

R2300BS

Plastic Polishing Cone Tool Base

DESCRIPTION: Tool base used to mount cone shaped sponge on.

R2260

Suction Cup (2 Piece)

DESCRIPTION: Two-piece suction cup used to hold onto contact lenses during modification. Has a convex and concave interchangeable piece to hold onto either side of the lens.

R2307

Round Sponge (3"x1" - Replacement)

DESCRIPTION: Round 3"x1" sponge with plexiglass base and screw set.

R2700

Drum Tool (3" - Complete)

DESCRIPTION: 3" velveteen covered tool with a flat 2 ¾" diameter sponge under the velveteen. Used for adding power, cleaning and polishing gas permeable and PMMA lenses.
INCLUDES: R2290 - 2 ¾" Edger sponge, R2700VP - 5" Velveteen pad, 3" Drum tool base

R2290

Edger Sponge (2 ¾")

DESCRIPTION: Flat 2 ¾" sponge used for 3" velveteen covered drum tool and also for automatic edgers (machines that automatically roll and polish edges on gas permeable and PMMA lenses.)

R2700VP

Velveteen Pads (5")

DESCRIPTION: Velveteen pad used for 3" drum tool.

R2630

Velveteen Pads (½") - (Not Pictured)