

# Bi expert<sup>®</sup>

GP Bifocal Design with




of comfort  
guaranteed.

**Bi**  
**expert**<sup>®</sup>  
Slab Off Technology

Surround your  
presbyopic patients in comfort  
with the new Bi expert GP bifocal design.

- Patented slab off technology – for a 360° uniform thin edge
- Non-truncated – resulting in a comfortable, thinner lens
- Rapid stabilization – delivering excellent visual clarity
- Available exclusively in these Boston materials – XO<sub>2</sub><sup>®</sup>, XO<sup>®</sup>, EO<sup>®</sup>, and ES<sup>®</sup>

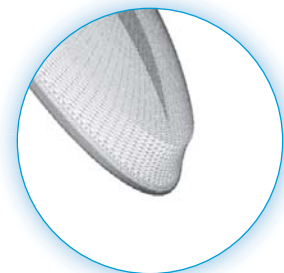

**Bi expert Design:**  
Inverse curve fits the lower  
eyelid for natural translation

*Our worry-free approach to successful GP fitting provides two refits  
per eye and full cancellation privileges for 120 days!*

Exclusively at

**ART**Optical  
contact lens, inc.

**BAUSCH + LOMB**  
Boston.

## GP Bifocal Design with Unique Edge Thickness Control


### Simple Lens Fitting – No Trial Set Needed

- Fit with K's and Rx, or can be more customized using other ocular measurements
- Immediate visual clarity for patients
- Less chair time to achieve final fit
- Fitting is easy with the Bi expert online calculator or the simplified Design Order Sheet – and the fit is guaranteed

### Key Design Features

- Unique slab-off construction naturally positions on lower eyelid\*
- Rapid, stable translation provides instantaneous vision change from distance to near
- Uniform edge thickness 360° around lens circumference, providing excellent comfort
- Available in spherical, aspherical, back toric, front toric designs
- Excellent option for all presbyopes

\*patent pending #PCT/FR2008/001125

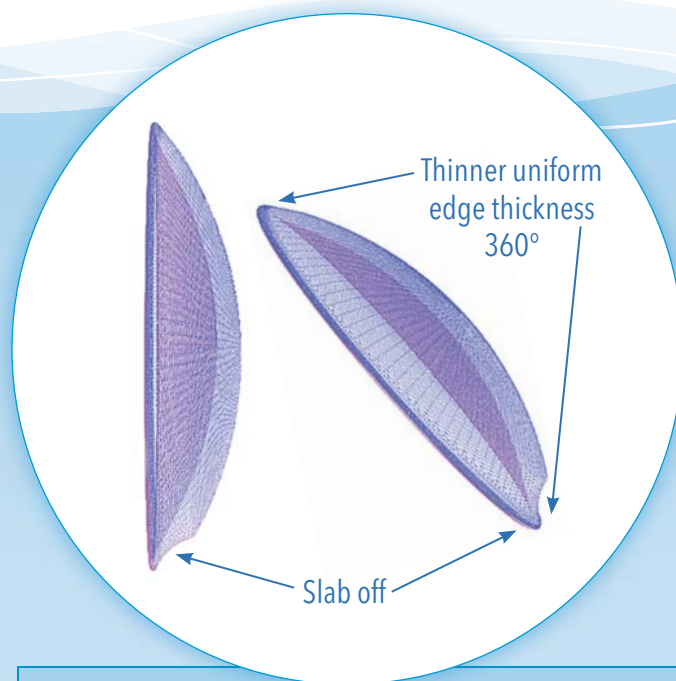

Lens Parameters	
Base Curve	7.20mm to 8.60mm
Diameter	8.50mm to 10.00mm
Distance Power	-15.00D to +10.00D in 0.25D steps
Add Power	+1.50D to +3.00D in 0.25D steps

Manufactured exclusively with Boston® Lens Materials

**BAUSCH+LOMB**  
Boston®



Dk 141



Dk 100



Dk 56



Dk 18

**ART**Optical  
contact lens, inc.

**Ordering** 800.253.9364

**Consultation** 800.566.8001

**Fax** 800.648.2272

**Online** [www.artoptical.com](http://www.artoptical.com)